DIETARIUS®
BE SMART BE HEALTHY

EN

USER MANUAL

Dietarius One D1120

We're delighted you've chosen the Dietarius One smart scale. Engineered with high-precision electrodes and advanced dual-frequency technology, this device provides you with comprehensive body-composition insights. Follow this guide to set up, operate, and maintain your scale for accurate, reliable measurements every time.

In our Dietarius app, you'll receive all of your essential measurement data at a glance—everything from weight and body-fat percentage to muscle mass and hydration levels is neatly organized and accessible. Simply sync your scale to the app to view detailed metrics, track trends over time, and gain personalized insights into your health journey.

TABLE OF CONTENTS

PACKAGING CONTENT	4
EQUIPMENT SYMBOL DESCRIPTION	5
SAFETY INFORMATION	6
THE MOST UP-TO-DATE VERSION OF THE USER MANUAL	9
PREPARATION FOR USE.....	10
INSTRUCTIONS	10
MEASUREMENT ON THE SCALE	11
SPECIFICATION PARAMETER INFORMATION	14
CLEANING & MAINTENANCE	16
TROUBLESHOOTING	17
BATTERY SAFETY AND ENVIRONMENTAL NOTICE	18
EU CONFORMITY	19
WARRANTY	20

- **Dietarius One smart scale**
- **USB-C cable**
- **User manual**

Trade name of manufacturer

Equitya s.r.o.

Address of manufacturer/importer

Silvanova 33, 90201 Pezinok ,Slovak Republic , EU

Product name

Dietarius-One

Batch, model, type

D1120

Product page

dietarius.eu/products/one

Country of origin

China

Production date

8/2025

Dietarius

is registered trademark of Equitya s.r.o.

DISCLAIMER INFORMATION IN THIS GUIDE MAY CHANGE WITHOUT NOTICE.

**Mobius loop / Recycling symbol**

the product or packaging is recyclable. please dispose of it responsibly at designated collection points.

**20 PAP**

paper material – used for cardboard packaging. recyclable in paper waste streams.

**04 LDPE**

low-density polyethylene – used in plastic bags or films. recyclable plastic.

**07 C/BE**

mixed or other plastics (e.g. bioplastics, multilayer materials). recyclability may vary; check local guidelines.

**Triangle with exclamation mark**

general warning – refers to important safety or cautionary instructions in the manual. always read carefully.

**CE**

this product complies with eu legislation for safety, health, and environmental protection requirements.

**Bluetooth® symbol**

indicates that the device includes bluetooth wireless technology for communication with smart devices.

**Wi-Fi symbol**

denotes presence of a wi-fi module, used for future features or connectivity (see manual for details).

**Lithium battery symbol**

the device contains a built-in lithium battery. do not open or dispose of in household waste. follow local battery disposal regulations.

**The crossed-out container symbol on the product**

in the accompanying documentation or on the packaging reminds you that in EU countries, all electrical and electronic equipment, batteries and accumulators must be disposed of separately at the end of their useful life as part of sorted waste. Do not dispose of these products in unsorted municipal waste.

Please read this manual carefully before use and keep it for future reference.

General Precautions

1. Read this manual carefully before first use and keep it for future reference.
2. This product is intended for household use only. Do not use it for medical, diagnostic, prevention of any disease (including obesity) or other professional purposes.
3. Children and persons with reduced physical, sensory, or mental capabilities may use the scale only under adult supervision.
4. Measurements are for informational purposes only and should not be taken as medical advice.
5. Ensure your feet are dry and clean before use.
6. Do not step onto the scale with wet feet or on wet surface of the scale, there is a risk of slipping.
7. For best consistency results, take measurements without clothing or footwear, prior to meals, and at the same time each day.
8. Do not use this product immediately after strenuous exercise.
9. Use the scale only for static pressure; do not jump onto or strike it.
10. Clean the scale with a damp cloth; do not use chemical solvents or immerse in water.
11. Place the scale on a flat, hard surface—never on carpet, uneven ground, or unstable platforms.
12. Do not place any objects on the scale when not in use.
13. This device must not be used during pregnancy.
14. Avoid standing on the edges of the scale to prevent risk of falling or injury.
15. Do not use the body fat analysis function if you have any of the following: a pacemaker, an electrocardiograph, or any other implanted or body-supported medical device.
16. Protect the scale from exposure to dust, chemicals, high humidity, toiletries, liquid cosmetic products, and abrupt temperature fluctuations. Avoid placing the device in direct sunlight.
17. Never disassemble, open, or expose the battery to moisture, high temperatures, or direct contact with conductive materials.
18. Only charge the device using the supplied cable and a certified adapter following the manufacturer's instructions

Prohibited Actions

1. Do not disassemble, crush, puncture, short-circuit, or expose the product to water, fire, or temperatures above 60 °C (140 °F).
2. Individuals with implanted medical devices (e.g., pacemakers) should not use this scale.
3. Avoid dusty, wet, dirty environments or strong electromagnetic fields to prevent internal circuit damage.
4. Operate only within 5 °C–40 °C (41 °F–104 °F) ambient temperature.
5. Keep the device away from heat sources or open flames.
6. Do not stand on wet or slippery surfaces.
7. Ensure your body and feet are dry before measuring to prevent slips, inaccurate readings, and water damage.
8. Do not attempt any internal repairs or modifications—only an authorized service center may perform maintenance.

Warnings

1. Never exceed 100 % of the rated weight capacity to avoid scale damage.
2. Stand in the center of the scale; do not stand on the edges or jump onto it.
3. Do not clean the product with hot water or volatile oils.
4. Keep away from radiators, stoves, fireplaces, etc., as it contains a lithium battery.
5. The product is not waterproof; avoid using in the rain or rinsing under running water.
6. Do not dismantle the battery — consult a repair center for all battery service.
7. The lithium battery is recyclable. Hand scale to an appropriate recycling facility when discarding the appliance.
8. The manufacturer shall not be held liable for any damage resulting from improper use of the device (e.g., damage to the platform or other components). Furthermore, the warranty is void if the safety precautions outlined above are not followed.
9. Please note that the data obtained from this device are indicative and may not match the precision of results derived from clinical medical examinations. Accurate assessment of biometric parameters can only be conducted by a qualified healthcare professional using certified medical diagnostic methods.
10. The Dietarius One smart scale contains a built-in 800 mAh lithium-ion battery, which includes a flammable electrolyte and reactive lithium compounds. If the battery is mechanically damaged, improperly charged, or experiences a short circuit, it may overheat, release toxic gases, or even ignite—potentially resulting in fire or explosion. This product complies with EU Battery Regulation (EU) 2023/1542 and Directive 2006/66/EC concerning lithium battery labelling, traceability, proper disposal, and safety. Please refer to the CE-marked label or manufacturer declarations for details.

Battery Care & Storage

- Use only the specified battery type and charger.
- If the scale will be unused for an extended period, charge the battery to at least 80% to prevent possible damage due to battery discharge and check the battery charge status at least once a month.
- Keep scale out of reach of children.
- Do not incinerate or dispose of batteries in household waste; follow local regulations for battery recycling.

Service & Maintenance

- No user-serviceable parts inside.
- All repairs, adjustments, or replacements must be carried out by an authorized service center.
- Failure to comply may void your warranty and create safety hazards.

Connection Requirement

11. To ensure proper functionality of the mobile application, an active internet connection is required. A stable Wi-Fi or mobile data connection (3G/4G/5G) is necessary for initial setup, account registration, synchronization of measurement data, and access to advanced features such as AI-powered insights and updates.
12. Please ensure your device is connected to the internet during use of the application to avoid interruptions or limited functionality.
13. Bluetooth Low Energy (BLE) connection is required to pair the scale with a compatible smart device (smartphone).
14. Prior to pairing, enable Bluetooth in your device settings and open the Dietarius app. The connection is secured via standard BLE pairing protocols
15. When you start scanning in the Dietarius app, it will automatically detect the Dietarius One scale once it's in pairing mode. Follow the on-screen instructions to complete the process.
16. All communication and data transfer via Bluetooth are fully automated. After launching the app and starting the scan, stand barefoot on the scale. The Bluetooth symbol will appear on the scale display, and data transmission to the Dietarius app will be triggered automatically.
17. For reliable connectivity, keep the scale and smart device within approximately 5 m (16 ft) during pairing and data sync.

Wireless Connectivity and Compatibility Notice

18. The manufacturer assumes no responsibility for any accidental, indirect, or consequential damages, including data loss or other harm resulting from information leakage during wireless communication via Bluetooth or Wi-Fi technologies.
19. Seamless wireless communication cannot be guaranteed with all smart devices. Devices must support Bluetooth specifications as defined by Bluetooth SIG, Inc., and must be appropriately certified for compliance. However, due to the wide range of devices available on the market, connectivity issues may still occur even with compliant devices. Such issues do not constitute grounds for a warranty claim.
20. Radio frequency transmissions from this device may interfere with the operation of medical equipment and other sensitive electronic systems. To avoid potential malfunctions or hazards, do not operate the device in the following environments:
 1. Near automatic doors, fire alarms, or other electronic safety systems
 2. Hospitals or other medical facilities
 3. Aircraft, trains, or other public transport systems
 4. Petrol stations or areas where flammable gases may be present

THE MOST UP-TO-DATE VERSION OF THE USER MANUAL

Due to DIETARIUS products and relevant regulations are regularly updated, the most up-to-date version of this user manual is always available online. Visit www.dietarius.eu or simply scan the QR code below using your smartphone to access the latest edition.

<https://dietarius.eu/products/one#manual>

PREPARATION FOR USE

Carefully remove all packaging materials and take the scale out of the box. Detach any protective films, stickers, or paper from the device. Charge the scale's battery using the supplied charging cable and a compatible charger. Once fully charged, place the scale on a flat, stable, and dry surface before use.

Proceed with the instructions in the App Installation & Setup section to complete the configuration.

INSTRUCTIONS

App Installation & Setup

1. Open the Apple App Store or Google Play Store and search for "Dietarius."
2. The Dietarius app may not be available in all regions.
3. Install the Dietarius app.
4. Make sure your device meets the minimum requirements before pairing:
 - o iOS: iPhone or iPad running iOS 15.0 or later
 - o Android: Version 9.0 or later with Bluetooth 4.0
5. Enable Bluetooth on your phone, launch the Dietarius app, and agree to the user terms.
6. Charge the scale fully before first use. Step on it gently to wake the display and confirm that Bluetooth (and GPS/location, if prompted) is active on your phone.
7. Sign in or register within the app, then enter your personal profile details.
8. In the app, go to BodyScan, then tap the green + icon in the top-right corner to add a new measurement.
9. Step onto the scale and wait for your phone to detect it—when the blue Bluetooth icon appears on the scale's screen, you're connected and ready to measure.

Working with the application is intuitive and it won't take you long to control it. Discover and enjoy the possibilities.

Note: This device is also equipped with a Wi-Fi module, the functionality of which is currently reserved for future features. These planned functions are not active at the time of product release and will be made available later through updates in the mobile application. The manufacturer reserves the right to modify the scope and availability of these features without prior notice.

Weighing Guidelines for Accurate Results

- Your scale is a precise electronic device. One of the most challenging objects to weigh accurately is the human body, as it is naturally subject to movement. For consistent and reliable measurements, always position your feet in the same place on the scale and remain as still as possible during weighing.
- For best results, weigh yourself at the same time each day—ideally in the morning—before eating or drinking. Keep in mind that heavy clothing, recent meals, or fluid intake can influence your weight. Short-term weight fluctuations are commonly due to changes in body fluid levels.

Basic Functionality Without the App

- When used without the DIETARIUS app on a smart device, the scale operates in basic mode and provides standard weight measurements, similar to a conventional scale.

Downloading the Application


- The app is available in both the Apple App Store and Google Play Store. To simplify the process, you may also scan the QR codes below to access the download page directly.
- Before using the device for the first time, please download the Dietarius application to your smart device free of charge. While the Dietarius app is free to use, certain advanced features are offered exclusively through paid subscription plans. An active subscription is required to access these premium services, which are optional and do not affect the basic functionality of the scale.

1. Preparation

- Remove all plastic film before first use.
- Always measure with clean, dry, bare feet in direct contact with the ITO electrodes.
- If you've just taken a reading, step off and wait for the display to show "0.00" before measuring again.


2. Proper Stance

Stand with your feet evenly spaced on either side of the scale's center. Your toes should rest on the front ITO electrode area and your heels on the rear half.

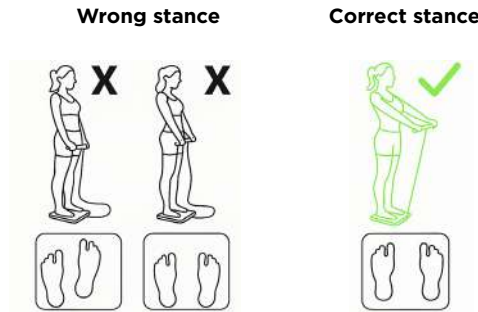3. Taking a Reading


1. Step onto the scale to begin recording.
2. When your weight stabilizes, the display will flash three times and then switch automatically to body-fat measurement.
3. Stay still until all values stop flashing and lock in place—they'll appear at the same time in the Dietarius app.
4. After about 10 seconds, the scale will shut off on its own.

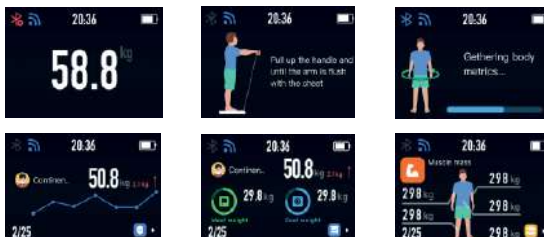Scale display during weighing.

4. Results Display

Upon completing a measurement, the Dietarius app automatically presents your results in a clear, organized layout. Each metric is grouped into intuitive categories, and interactive charts let you visualize your progress and trends over time. To help you make sense of the data, the Dietarius AI Coach analyzes your results and provides personalized, actionable insights.

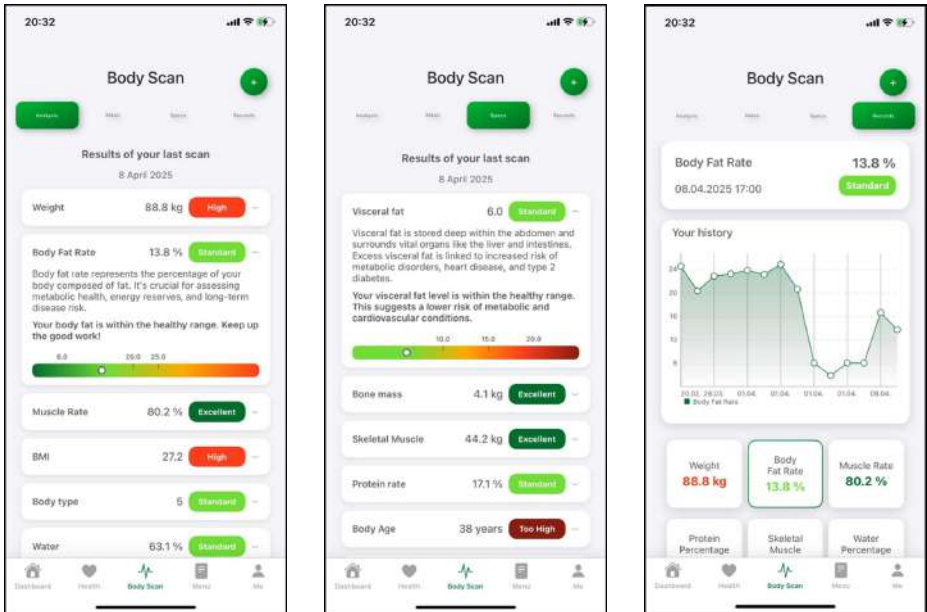Measurement result overview in Dietarius app

Precautions

- The product is activated by pressure (by pressing or standing on the scale)
- After it is placed on a hard and flat ground, it will turn on the screen after sensing a weight of more than 3 kg.
- For products with built-in lithium batteries, please use a 5 V-1 A charger to charge for 2 hours when using for the first time.
- After using the scale, put it into the horizontal position (not vertically), out of reach of children and incapacitated people (see par. I. SAFETY INFORMATION).

The Dietarius app processes the following items

Weight · Body Fat · BMI · Muscle Mass · BMR · Recommended Calorie Intake · Body Water · Fat Mass · Protein · Moisture Content · Protein Amount · Lean Body Mass · Subcutaneous Fat Rate · Skeletal Muscle Rate · Visceral Fat · Muscle Ratio · Muscle Control Quantity · Fat Control Quantity · Standard Weight · Weight Control · Subcutaneous Fat Mass · Left-Hand Fat Percentage · Body Cell Mass · Right-Hand Fat Percentage · Intracellular Water Mass · Trunk Fat Percentage · Extracellular Water Mass · Left-Leg Fat Percentage · Inorganic Salt Content · Right-Leg Fat Percentage · Left-Hand Muscle Mass · Left-Hand Fat Mass · Right-Hand Muscle Mass · Right-Hand Fat Mass · Trunk Muscle Mass · Trunk Fat Mass · Left-Leg Muscle Mass · Left-Leg Fat Mass · Right-Leg Muscle Mass · Right-Leg Fat Mass · Obesity Rate · Bone Mass · Body Type · Obesity Level · Health Assessment · Body Age · Body Score

SMART Device Dietarius One

Intelligent Body-Fat Scale with 8 Advanced-Sensor Electrodes

- **Product Name:** Dietarius One
- **Model Number:** D1120
- **Product Dimensions:** Product Dimensions (L × W × H): 335 × 335 × 33 mm
- **Measurement Range:** min. 3 kg - max. 180 kg
- **Capacity:** max.180 kg / 400 lb
- **Resolution:** 0.05 kg / 0.1 lb
- **Platform Material:** 6 mm Tempered Glass
- **Display:** 3.5-inch Color TFT
- **Charging Port:** USB Type-C
- **Scanning time:** up to 30 seconds
- **Battery:** 800 mAh lithium batteries
- **Certifications:** CE · RoHS · FCC

Connectivity & Module

- Low-energy Bluetooth module, compatible with iOS and Android
- Real-time dual-mode updates via Bluetooth and (future) Wi-Fi

Wireless information

Mode	Frequency Band (MHz)	Maximum Output Conducted and Power (dBm)
BT LE	2402-2480	5 dBm
WLAN	2412-2462	16 dBm

Key Features

- App processes up to 56 parameters (see App. processes the following items) : weight, standard weight, body fat, body water, BMI, BMR, bone mass, muscle mass, visceral fat, protein, fat rate, body type, subcutaneous fat rate, lean body mass, etc.
- 8 High-Precision Sensors for accurate pressure and impedance readings
- Step-on Activation & Auto-Power-Off to conserve battery
- Low-Battery & Overload Indication

Advanced Functionality

The Dietarius One dual-frequency scale embeds a advanced analytical algorithm that automatically calculates and analyzes body-composition metrics from its eight electrodes. In addition to total body measures, segmented readings for arms, legs, and trunk provide a more detailed regional analysis.

- Heart-Rate Test on Screen: integrates a quick pulse check for cardiovascular insight alongside body composition.
- VA-Type Display: delivers crisp, detailed visuals of all collected data, helping users track progress at a glance.
- AI-Based Measurement Chip: ensures high measurement precision and stable readings
- ITO-Coated Electrodes: standard on body-fat scales to guarantee consistent conductivity and precise impedance measurement.

Power & Endurance

An 800 mAh rechargeable Li-ion battery powers the device for extended use. After a brief 2-hour initial charge via the Type-C port, the scale is ready for continuous operation. Maximum supported user weight is 180 kg.

Note: Bluetooth connectivity is fully operational. Wi-Fi support will be available in a future firmware and app. update.

Power Supply Requirements

To ensure safe and reliable charging, use a power adapter that provides a voltage of 5 V and a current of at least 1 A, corresponding to a minimum power output of 5 W. This meets the device's electrical requirements.

It is recommended to use a certified USB-C charger that complies with international safety standards (e.g. CE, UL). Using low-quality or uncertified chargers may result in improper charging, overheating, or potential damage to the built-in lithium battery.

CLEANING & MAINTENANCE

To keep your scale in optimal condition, gently wipe the surface with a soft sponge moistened with warm water, then dry thoroughly with a clean cloth.

Warning:

- Never immerse the scale in water or any other liquid.
- Do not use abrasive pads (plastic or metal), harsh detergents, chemicals, gasoline, solvents, or similar substances, as these can severely damage the finish.

Storage

Before storing, ensure the scale is thoroughly clean and dry. Store it in a well-ventilated, dry location out of reach of children and pets. Do not place any items on top of the scale during storage.

If you experience any issues, try the solutions below. If problems persist, please contact your supplier.

- Scale won't power on – Ensure the scale is charged or properly connected to its power source.
- Display shows “FAIL” – The Wi-Fi setup failed. Open the app and reconfigure your network settings.
- Inaccurate or no data during measurement – Make sure you're standing with clean, dry, bare feet directly on the electrodes. – If body-fat readings fail, step off and try again.
- “!” Overload warning – You have exceeded the 180 kg (400 lb) capacity. Please step off immediately.
- Blank screen, errors, or missing icons – Verify that Bluetooth is enabled on your phone.
- Sudden shutdown during use – The scale has built-in overheating protection. If it turns off, let it cool for a few minutes before restarting.

Battery Safety and Environmental Notice

Warning: This product contains a built-in 800 mAh lithium battery. For your safety, do not attempt to remove, replace, or repair the battery yourself under any circumstances. Improper handling of lithium batteries may lead to fire, explosion, or injury.

Battery Information: The internal lithium battery is designed for long-term use and does not require user maintenance. Charging should only be performed using the supplied USB-C cable and a compatible certified power adapter.

Environmental Responsibility and Disposal: This product contains electronic components and a lithium battery and must not be disposed of with household waste. At the end of its life, **please return the device to a designated collection point for the recycling of electrical and electronic equipment**, in accordance with local regulations (e.g. WEEE Directive in the EU).

By disposing of this product responsibly, you contribute to environmental protection and the sustainable use of natural resources.

To support sustainability and responsible consumption, this product and its packaging may include labels indicating the types of materials used and their recyclability, where space permits.

Symbols displayed on the product or in its accompanying documentation indicate that electrical and electronic equipment must not be disposed of with regular household waste. Instead, please take such items to designated electronic waste collection points, where they will be accepted free of charge. Proper disposal helps preserve valuable natural resources and prevents potential harm to the environment and human health that may result from inappropriate handling of electronic waste.

Battery Disposal Notice

The Dietarius One smart scale contains a built-in lithium battery, which includes substances that may be harmful to the environment or health. Do not dispose of this device or its battery with regular household waste.

Please return the used device to an authorized collection or recycling facility for proper handling and disposal in accordance with local environmental regulations.

For detailed information on recycling and disposal, contact your local waste management authority or collection center. Be aware that improper disposal may result in legal penalties under national regulations.

Servicing and Warranty:

Any maintenance that involves internal components must be performed by authorized service professionals only. Failure to comply with the manufacturer's instructions will void the warranty or lost guaranty for repair

EU conformity:

CE Declaration of conformity

Equityta s.r.o. hereby declares that the Dietarius One product meets the requirements of the standards and regulations that are relevant for the type of equipment, including the RED Directive (2014/53/EU).

Equityta s.r.o. hereby declares that the device Dietarius One Smart scale complies with the relevant EU legislation and international requirements in force, including: - Battery Regulation (EU) 2023/1542 for internal lithium batteries; - Radio Equipment Directive 2014/53/EU and EMC Directive 2014/30/EU for wireless connectivity (Bluetooth and Wi-Fi); - RoHS Directive 2011/65/EU restricting hazardous substances; - WEEE Directive on responsible recycling and disposal. The product bears the CE marking in accordance with applicable Union harmonisation legislation.

A full EU Declaration of Conformity including applied harmonised standards is available on <https://dietarius.eu/products/one/>

Warranty

Dietarius® Two (2) Years Limited Warranty

Equitya s.r.o., Silvanova 33, 902 01 Pezinok, Slovak Republic, EU (“Equitya”) warrants the Dietarius One-branded hardware product (“Dietarius Product”) against defects in materials and workmanship when used normally and in accordance with Equitya’s published guidelines for a period of two (2) years from the date of original retail purchase by the end-user (“Warranty Period”).

Equitya’s published guidelines include, but are not limited to, information contained in technical specifications, safety instructions, or in this guide.

This warranty does not guarantee uninterrupted or error-free operation of the Dietarius Product. Equitya is not responsible for damage arising from failure to follow instructions relating to the proper use of the Dietarius Product.

Manual version 1.3 08/2025

Made in China

Copyright © Equitya s.r.o.

www.dietarius.eu

All rights reserved.

The appearance and technical specifications of the product are subject to change without notice.

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